


# CUTISCOPE USER MANUAL

Model No: DRM100

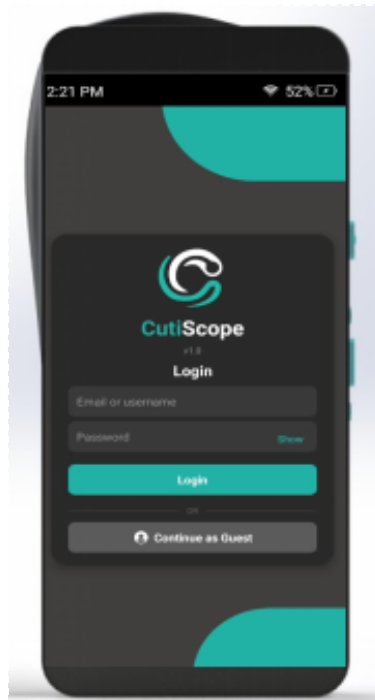

**Read all instructions and safety warnings before use. Ensure the device is fully charged and disinfected before its first clinical application.**

## Table of Contents

|                                                               |           |
|---------------------------------------------------------------|-----------|
| <b>1. Introduction.....</b>                                   | <b>1</b>  |
| 1.1. Intended Use.....                                        | 1         |
| 1.2 User Profile.....                                         | 1         |
| 1.3 Warnings.....                                             | 1         |
| 1.4 Safety Precautions.....                                   | 2         |
| <b>2. Technical Specifications.....</b>                       | <b>3</b>  |
| <b>3. Device Components.....</b>                              | <b>4</b>  |
| 3.1. Device Structure.....                                    | 4         |
| <b>4. Device Overview.....</b>                                | <b>5</b>  |
| 4.1. Device Operating Instructions.....                       | 5         |
| 4.2 Device Charging and Disinfection.....                     | 6         |
| <b>5. System Overview.....</b>                                | <b>7</b>  |
| <b>5.1. Authentication &amp; Security.....</b>                | <b>8</b>  |
| 5.1.1. Clinician Sign-In.....                                 | 8         |
| 5.1.2. Guest Session Protocol.....                            | 9         |
| <b>5.2. Camera Screen(Clinician Interface).....</b>           | <b>10</b> |
| 5.2.1. Patient Addition.....                                  | 11        |
| 5.2.2. Image Capturing.....                                   | 12        |
| 5.2.3. Gallery.....                                           | 15        |
| 5.2.4. User settings.....                                     | 19        |
| <b>5.3. Wi-Fi Connectivity.....</b>                           | <b>20</b> |
| <b>5.4. System Calibration (Date &amp; Time).....</b>         | <b>21</b> |
| <b>5.5. Power ON/OFF.....</b>                                 | <b>22</b> |
| <b>5.6. Sleep Mode.....</b>                                   | <b>23</b> |
| <b>5.7. Battery Charging.....</b>                             | <b>23</b> |
| <b>5.9 Web Dashboard.....</b>                                 | <b>24</b> |
| 5.9.1. Sign Up /Log In.....                                   | 24        |
| 5.9.2. Dashboard.....                                         | 25        |
| 5.9.2.1. Image Gallery.....                                   | 25        |
| 5.9.2.2. Patients.....                                        | 26        |
| 5.9.2.3. Audit Logging.....                                   | 26        |
| <b>6. Electrical Safety.....</b>                              | <b>27</b> |
| <b>7. Cleaning Protocol (Routine - Before/After Use).....</b> | <b>28</b> |
| <b>8. Troubleshooting Guide.....</b>                          | <b>29</b> |
| <b>9. Manufacturer Details.....</b>                           | <b>31</b> |

## 1. Introduction

The Revive Meditech **Cutiscope** is a handheld contact-based hybrid dermatoscope designed to provide enhanced visualization and analysis of skin structures using both polarized and non-polarized illumination.

The system includes:

- Dermatoscope device
- Simple push-button controls
- Rechargeable battery with USB-C charging
- Image capture software
- Integrated optical system (polarized & non-polarized modes)
- CutiScope Cloud

With up to **10x optical and 3x digital magnification**, the Cutiscope allows precise observation of skin features.

### 1.1. Intended Use

The device is intended for non-invasive examination of skin lesions and image capturing to assist healthcare professionals in dermatological evaluation.

### 1.2 User Profile

This device is intended for use by:

- Dermatologists
- Physicians
- Trained healthcare professionals

### 1.3 **Warnings:**

- Do not use the device for diagnosis without proper clinical training or supervision.
- Do not use the device with wet hands.
- Do not use the device on open wounds, infected skin, or highly sensitive areas unless advised by a qualified healthcare professional.
- Avoid direct exposure of the device light to the eyes, as it may cause discomfort or injury.
- Do not attempt to disassemble or modify the device; this may result in malfunction or electrical hazards.

- Use only the provided or recommended charging accessories to prevent damage or risk of fire.
- Do not operate the device in wet or highly humid environments to avoid electrical hazards.
- Ensure the device is properly cleaned before and after each use to prevent cross-contamination.
- Do not operate the device if the device overheats, emits an unusual odor, or shows signs of damage.
- Keep the device out of reach of children and untrained individuals.

#### 1.4 Safety Precautions:

| Category                 | Precaution                                                                                                                                                                                          |
|--------------------------|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| <b>General Use</b>       | Operate the device only as per the user manual instructions.                                                                                                                                        |
| <b>Handling</b>          | Handle the device carefully, avoid mechanical shock or drops as it contains sensitive optical and electronic components.                                                                            |
| <b>Hygiene</b>           | Clean and disinfect the contact surface before and after each use. Use NON ALCOHOLIC WIPES to disinfect. Do not use IPA or Alcohol-based disinfectants.                                             |
| <b>Charging</b>          | <ol style="list-style-type: none"> <li>1. Use only approved USB-C chargers and avoid overcharging.</li> <li>2. Do not operate the device while it is charging.</li> </ol>                           |
| <b>Environment</b>       | Use the device in a clean, dry, and well-lit environment.                                                                                                                                           |
| <b>Storage</b>           | Store the device in a protective case when not in use.                                                                                                                                              |
| <b>Temperature</b>       | Avoid exposure to extreme heat or cold conditions.                                                                                                                                                  |
| <b>Battery Safety</b>    | Do not puncture, crush, or expose the battery to high temperatures.                                                                                                                                 |
| <b>Light Safety</b>      | <ol style="list-style-type: none"> <li>1. Avoid exposure to light in the eyes. Do not use it for eye examinations.</li> <li>2. Ensure patients keep their eyes closed during facial use.</li> </ol> |
| <b>Infection Control</b> | Follow proper disinfectant protocols to prevent cross-patient contamination.                                                                                                                        |

## 2. Technical Specifications

| Item                                    | Specification                                                                                                  |
|-----------------------------------------|----------------------------------------------------------------------------------------------------------------|
| <b>Net Weight/ Dimension</b>            | 348.5 g/ 175x84x45 mm                                                                                          |
| <b>Package Weight/ Dimensions</b>       | 607.5 g/ 200x105x67 mm                                                                                         |
| <b>Operating and Storage Conditions</b> | 10°C to 40°C, 30%–75% RH (non-condensing)                                                                      |
| <b>Light Source</b>                     | LEDs                                                                                                           |
| <b>LED Type</b>                         | SMD LED Beads                                                                                                  |
| <b>Polarization</b>                     | 15 LEDs (10 polarizing and 5 non-polarizing)                                                                   |
| <b>Aperture</b>                         | 45 mm                                                                                                          |
| <b>Digital Magnification</b>            | 3x                                                                                                             |
| <b>Optical Magnification</b>            | 10x                                                                                                            |
| <b>Function buttons</b>                 | <ul style="list-style-type: none"> <li>• Mode switching button.</li> <li>• Power ON/OFF Push button</li> </ul> |
| <b>Examination Method</b>               | Contact only                                                                                                   |
| <b>Battery</b>                          | Rechargeable Lithium-Ion Battery                                                                               |
| <b>Charging Time</b>                    | Min. 3 hours                                                                                                   |
| <b>Continuous duty time</b>             | Approx. 2.5-3 hours                                                                                            |
| <b>Charging Terminal</b>                | USB Type C                                                                                                     |
| <b>Wifi Support</b>                     | Yes                                                                                                            |


### 3. Device Components

Ensure the following items are included:


- Dermatoscope Unit
- USB Type-C charging cable
- Power adapter
- Electronic user manuals can be accessed by QR Code on devices.

If any component is missing, contact the manufacturer (refer to [Section 9: Manufacturer Details](#)).

#### 3.1. Device Structure:

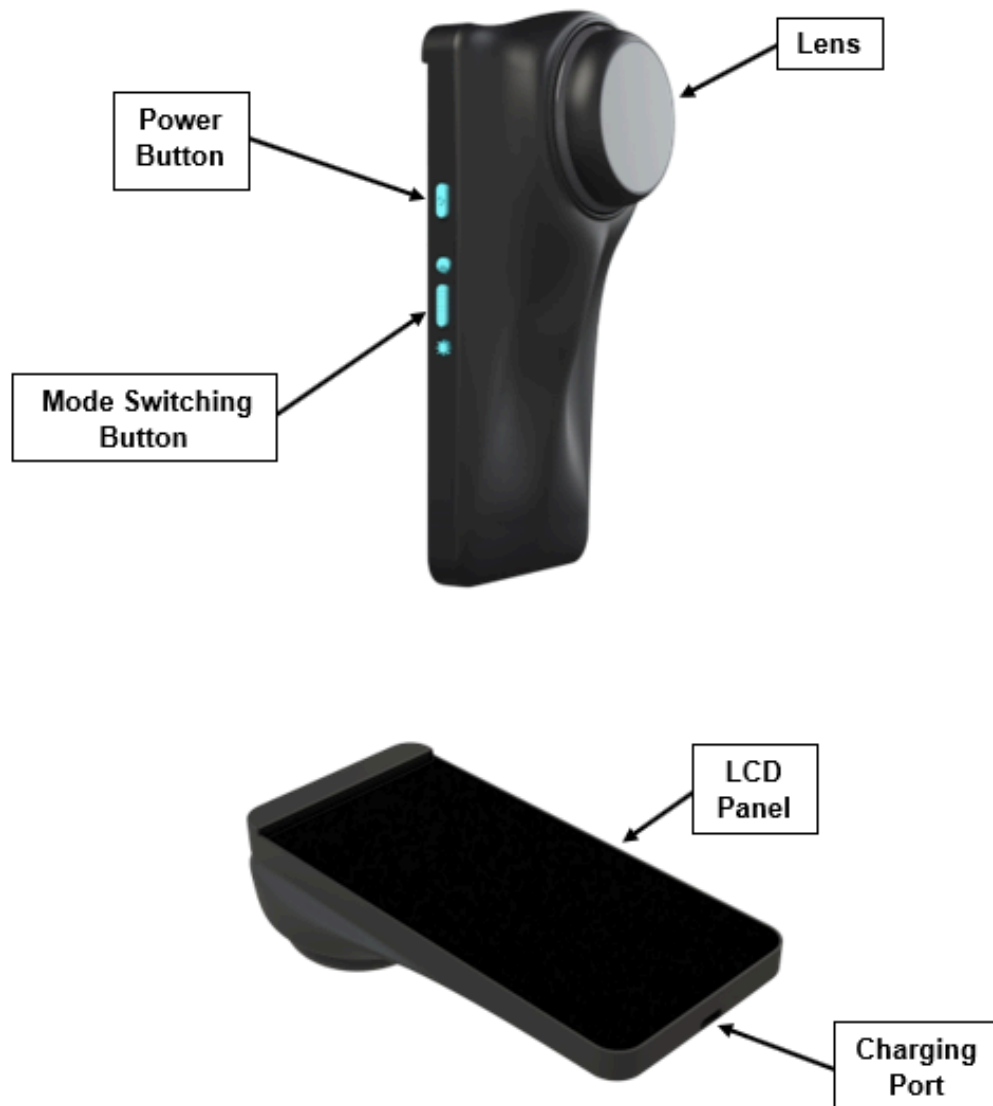

## 4. Device Overview

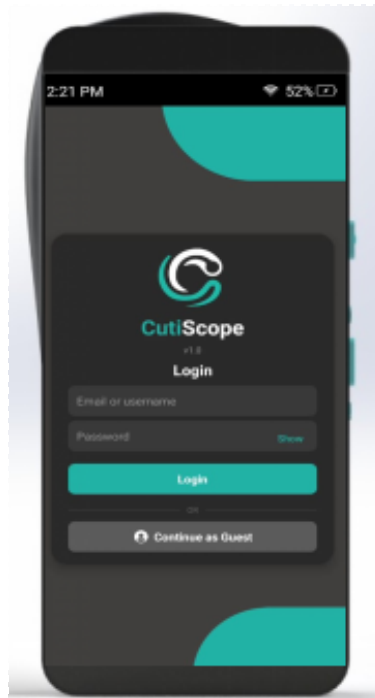

### 4.1. Device Operating Instructions

The device should be used according to proper operating guidelines to ensure safe and efficient performance. Follow the steps below:

#### 1. Charging the Device

- Connect the charger to the **charging port** on the base unit.
- Ensure the device is fully charged before use.

#### 2. Powering On

- Press and hold the **power button** to turn ON the device.
- The **LCD panel** will light up, indicating the device is ready.

#### 3. Illumination & Mode Activation

- Illuminate the LEDs using the illumination switch present on the camera screen. 
- Use the **mode switching button** to toggle between:
  - **Polarized Mode** 
  - **Non-Polarized Mode** 
- The selected mode will be indicated through the **mode indicator on the dermoscope application** (refer to [Section 5.2. Camera Screen](#)).

#### 4. Operation

- Position the **lens** over the target area.

- Switch on the illumination in the app.
- Adjust modes if required.
- Observe and capture images using the application.

#### 5. Device Handling & Usage Limitation

- Avoid operating the device continuously for periods exceeding **40 minutes**.
- After each usage cycle, allow the device to **cool down** before reusing.

#### 6. Powering Off

- Press and hold the **power button** to turn OFF the device.

#### 7. After Use

- Recharge the device if needed.
- Store in a safe, dry place to protect the **lens** and components.

## 4.2 Device Charging and Disinfection

Proper charging and cleaning of the device is necessary for long-term durability and hygiene.


- Charge the device using the recommended charger only
- Avoid overcharging
- Keep the charging port clean and dust-free
- Disinfect the device regularly using a soft cloth and Alcohol-free disinfectant
- Avoid excessive liquid contact with the device.



## 5. System Overview:

The Cutiscope Device operates on **software** designed to ensure reliability and focus during examinations connected with a cloud storage **web application**.


- **Unified Operational Environment:** The system is dedicated exclusively to CutiScope Dermatoscope image capture software.
- **Auto-Standby:** To preserve power and clinician privacy, the CutiScope Dermatoscope enters a standby state after 2 minutes of inactivity.

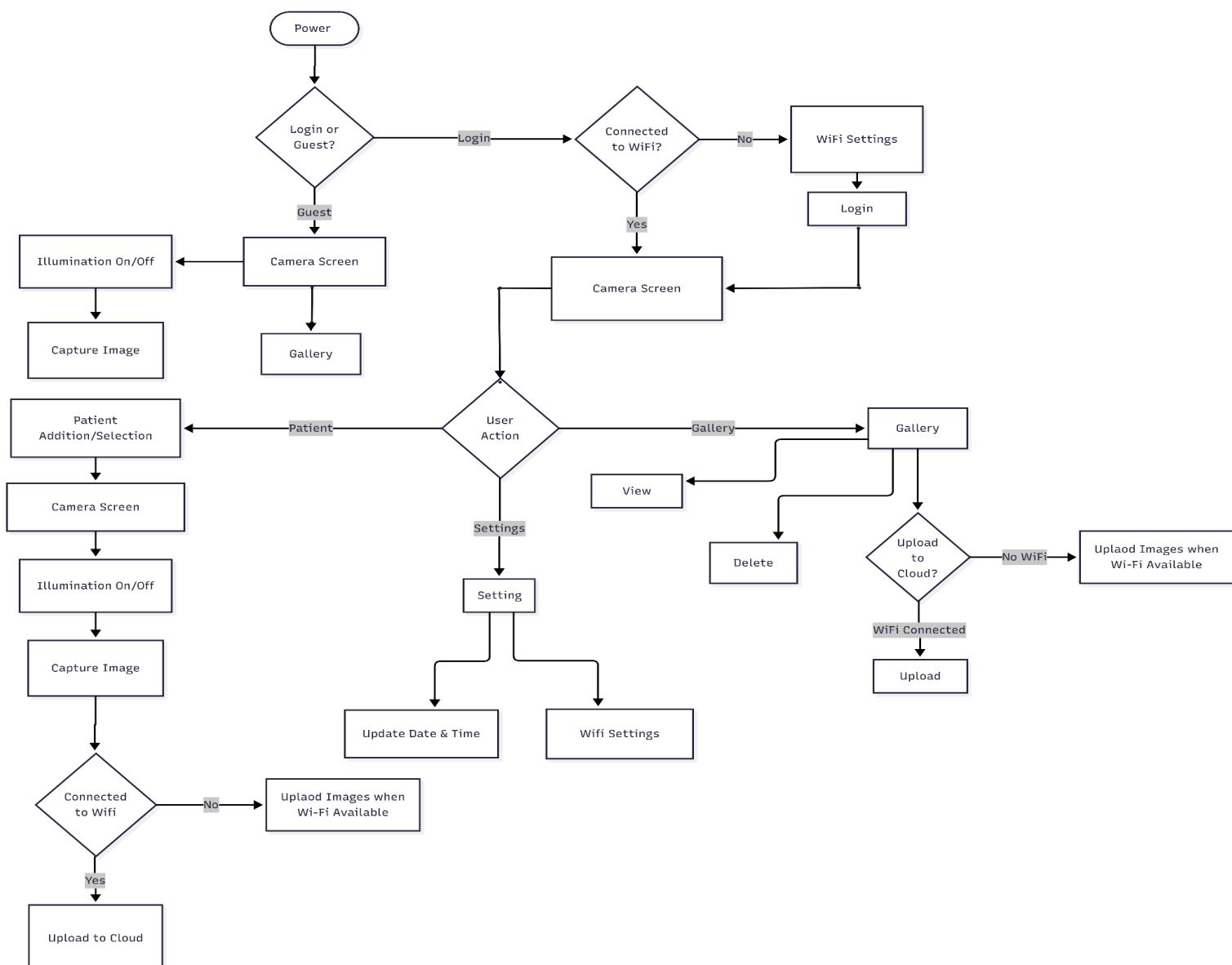


## 5.1. Authentication & Security

### 5.1.1. Clinician Sign-In:

- Authorized medical professionals must authenticate their credentials. This ensures every diagnostic capture is linked to a specific clinician for accurate record-keeping and auditing.
- The login credentials that are used will be the same that clinician use to sign-up in the web dashboard (refer to [Section 5.9.1: Sign Up / Log In](#) ).

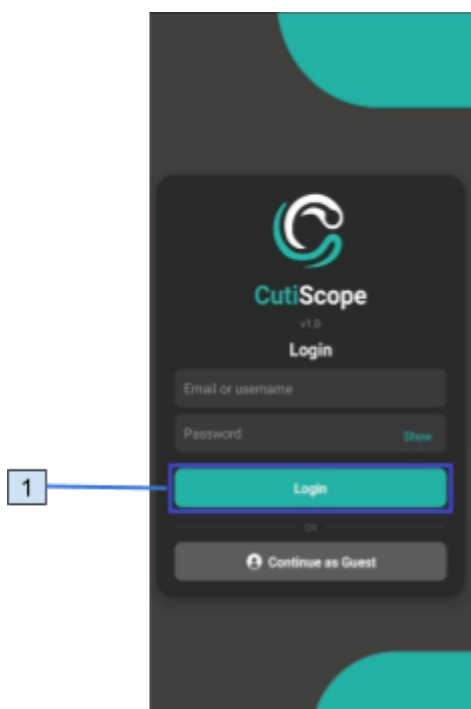

Login as Clinician

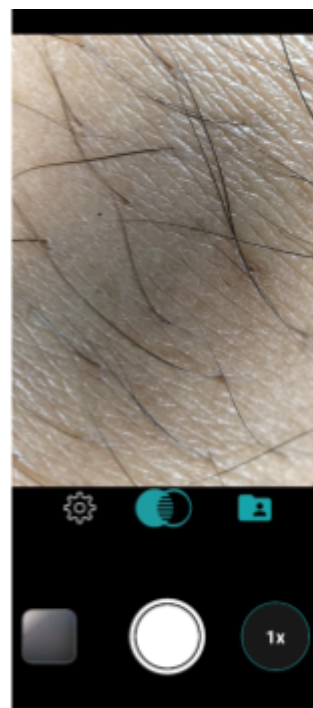


Camera Screen


1.Login button

### 5.1.2. Guest Session Protocol:


A temporary Guest session is available for immediate use without a permanent account. The images that are captured during the guest mode will be saved temporarily and will be removed upon exit.

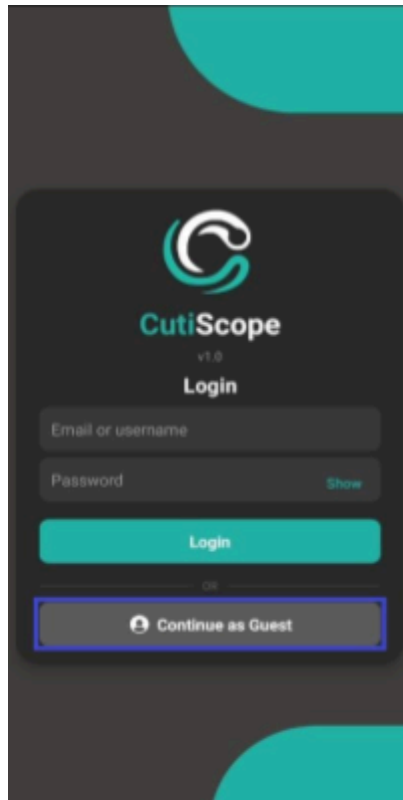

Continue as Guest

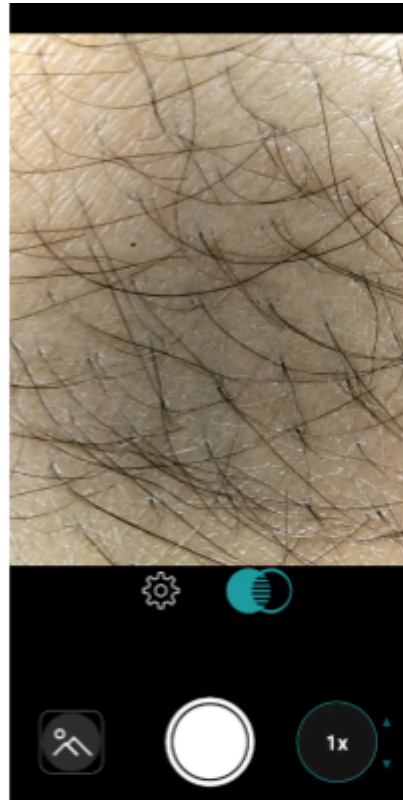

Camera Screen

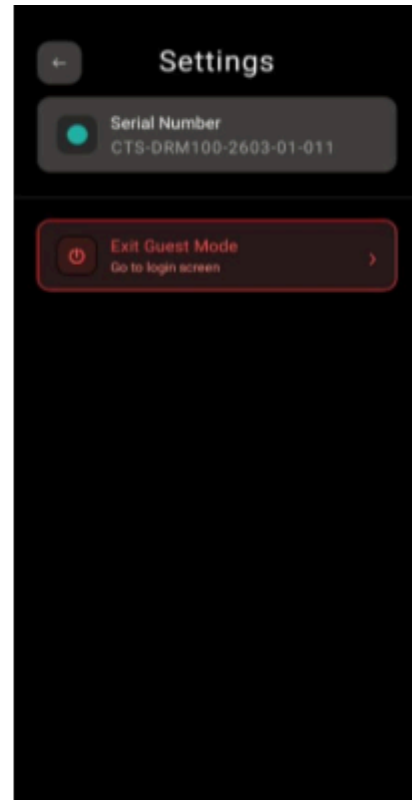

Settings Screen


### ⚠ CAUTION:

**SESSION DATA PERMANENCE:**For patient privacy, all data captured during a Guest session is stored in volatile memory and is permanently erased upon session termination.

## 5.2. Camera Screen(Clinician Interface):

In the camera screen, the user can navigate to other screens, such as:

1. Patient Addition
2. Image Capturing
3. Gallery
4. Settings

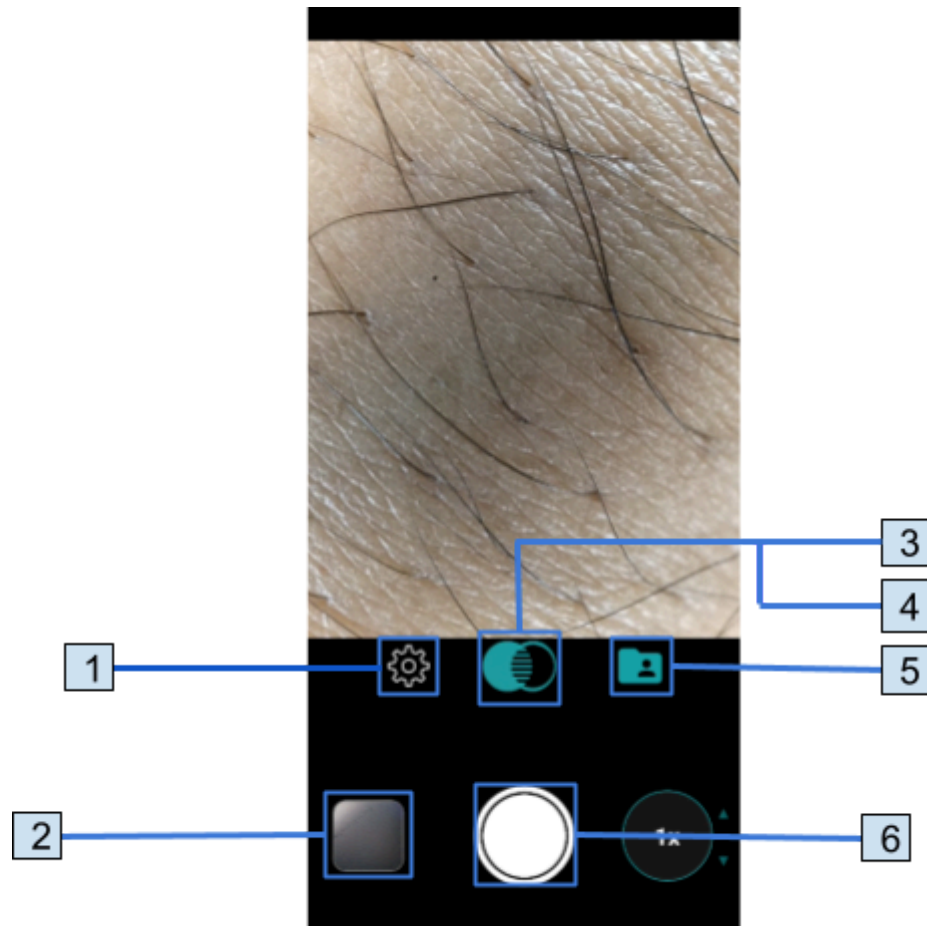


- 1.Settings
- 2.Gallery
- 3.Mode Indicator
- 4.Illumination Switch
- 5.Patient Addition
6. Capture Image

-  Off State
-  Polarization Mode
-  Non-Polarization Mode

### 5.2.1. Patient Addition:


- When a clinician logs in, first he/she will add a new patient or select an existing one to capture an image.
- The captured image would be saved with the selected patient's name and in the folder created for that patient in the gallery to store captured images.
- The patients that are present in the list would be added from the web dashboard (refer to [Section 5.9.2.2 : Patients](#))

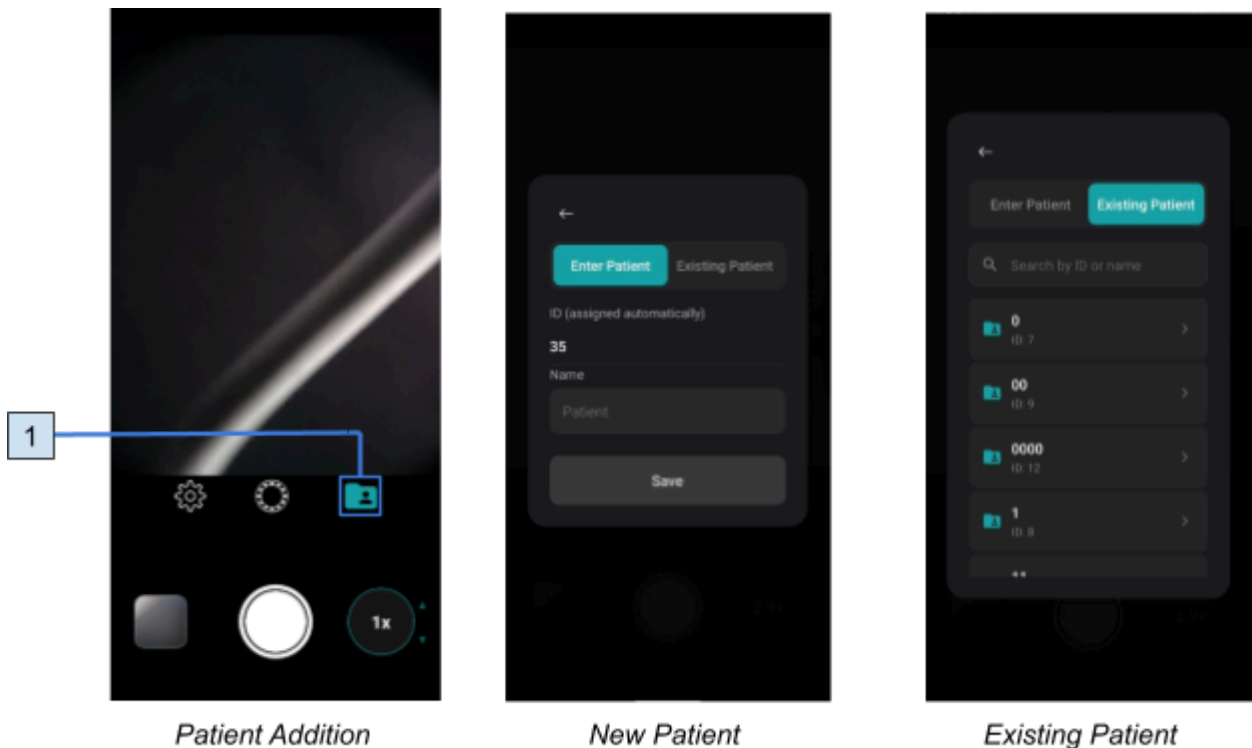

#### 1.Patient Addition

### 5.2.2. Image Capturing:

- The image will be captured only if the user selects the patient from the list or enters a new patient as the captured images will be saved to that particular folder.
- Pressing the center button activates the Lens LED in polarized mode, accompanied by an on-screen indication.
- To toggle between polarized and non-polarized modes, press the lower button on the left side. The indication will update accordingly to reflect the selected mode.

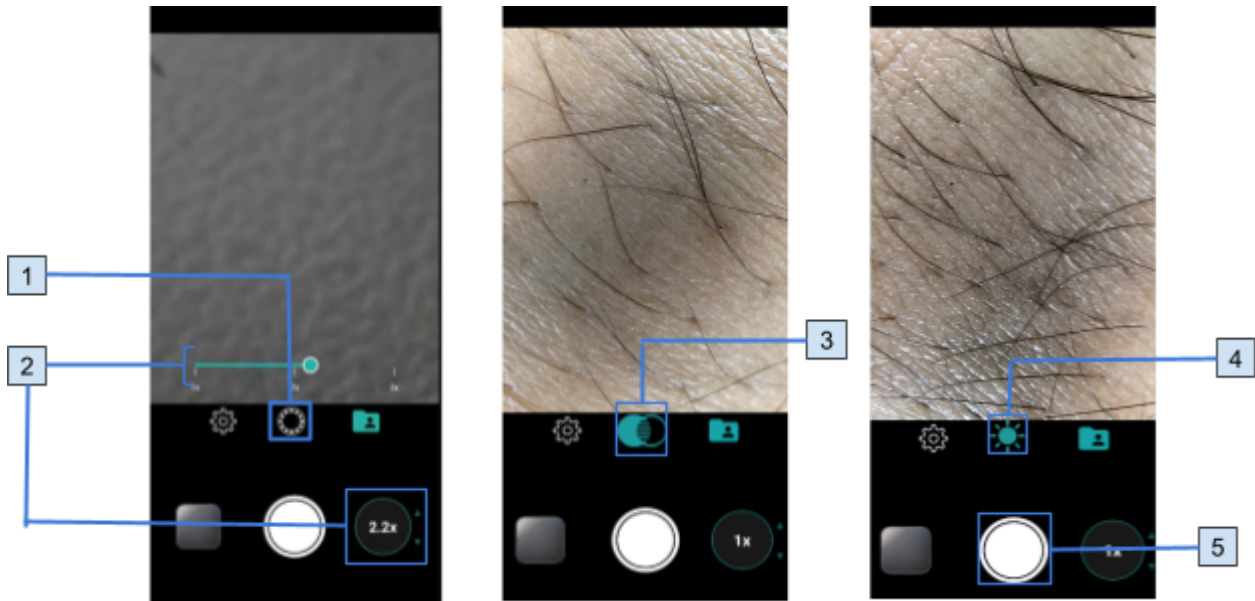


1. Off State
2. Zoom Control
3. Polarization
4. Non-Polarization
5. Capture Image

- As Cutiscope is a contact-based dermatoscope, place the device gently on the skin, ensuring proper contact.
- Adjust focus and magnification.


**Zooming:**

- Use the slider to fine control the zooming.
- Adjust magnification using on-screen gestures
- Use the slider to fine control the precision adjustment.

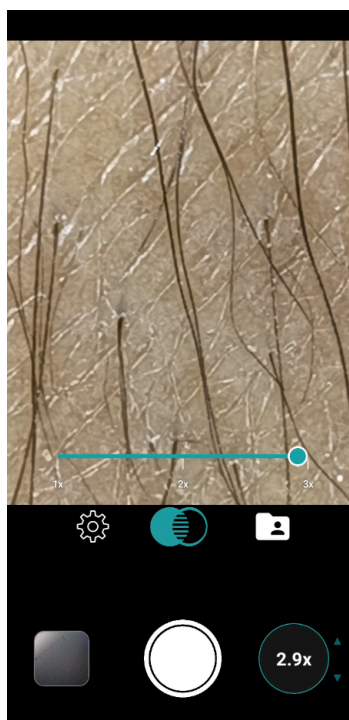

**Focus:**




- The lens has both auto focus and manual focus features.
- Tap the screen to focus manually.
- Capture images using a software interface.

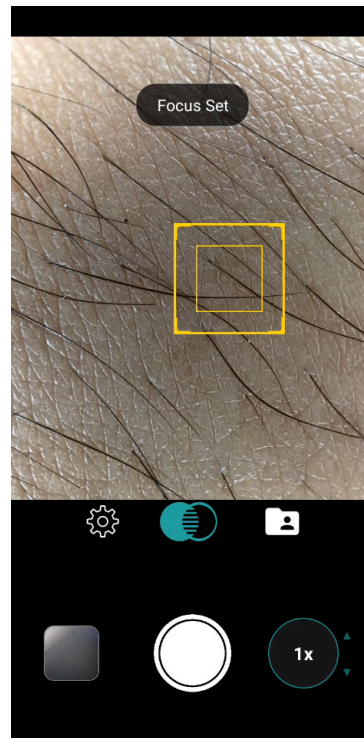







### 5.2.3. Gallery:

- In the gallery, the images that are captured will be saved in the folders.
- A **main folder** is created using the **Patient Name**.
- Inside the patient folder, **sub-folders are created for each Year** (e.g., 2026, 2027).
- Within each year folder, **monthly folders** are created (e.g., January, February).
- All captured or uploaded **images are stored inside the corresponding month folder**.

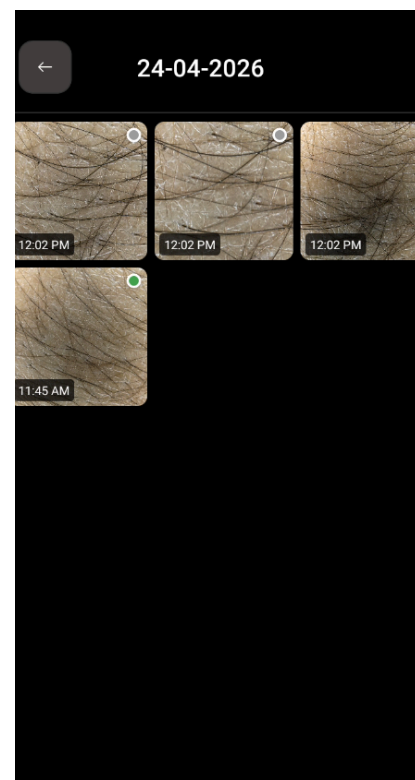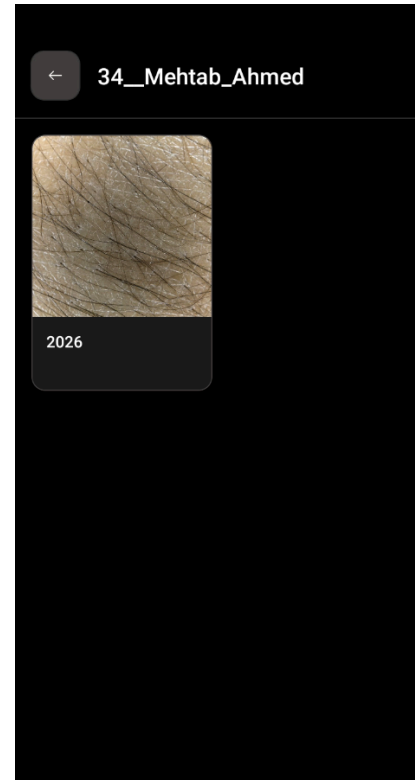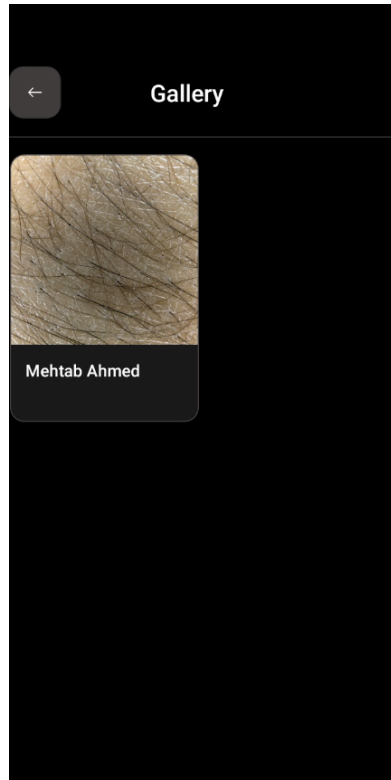




### Image Uploading:

- Images that are captured are uploaded to the cloud if the Internet is available(refer to [Section 5.9.2.1 : Image Gallery](#)).
- If the Internet is not available, then the images are stored in the gallery, but the user will have to upload the image manually, either by selecting all at once or one at a time.
- The stored images are seen on the web app in the specific folder of the patient's name(refer to [Section 5.9.2.1 : Image Gallery](#)).
- The images that have a green indicator on the top-right corner are already uploaded to the cloud, and the remaining are not.

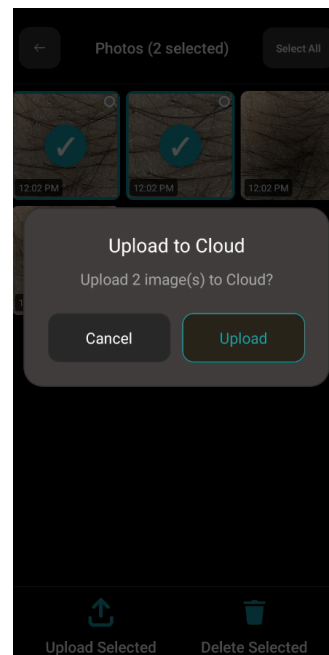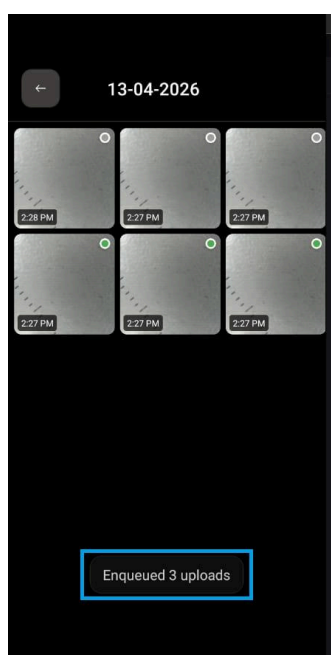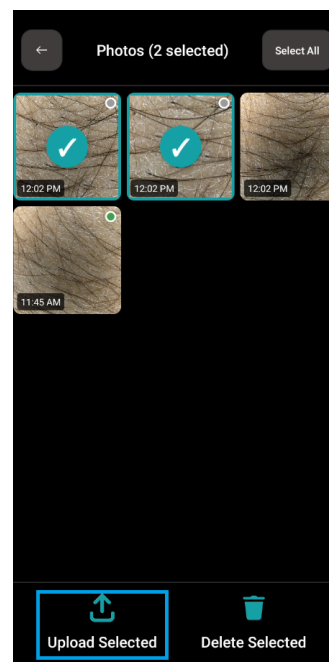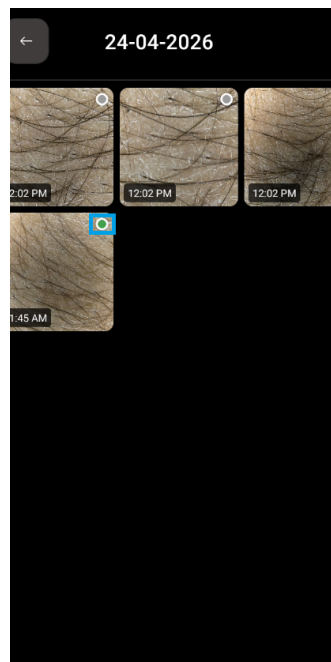

### Image Deletion:

- Images that are captured will be deleted from the gallery all at once by selecting the album.
- The user can delete the images by selecting one by one or all at once.
- The images that are uploaded to the cloud will be deleted from the device but not from the cloud storage.

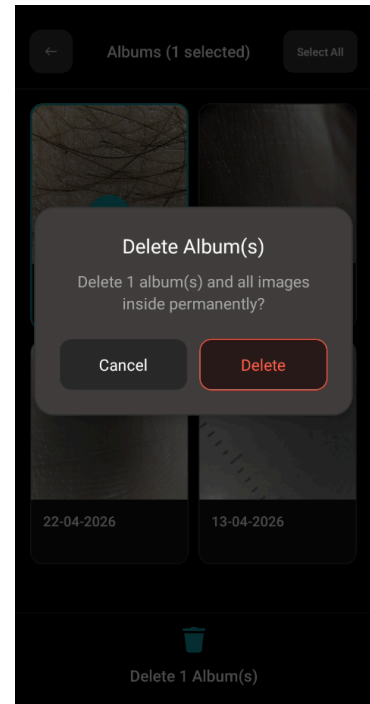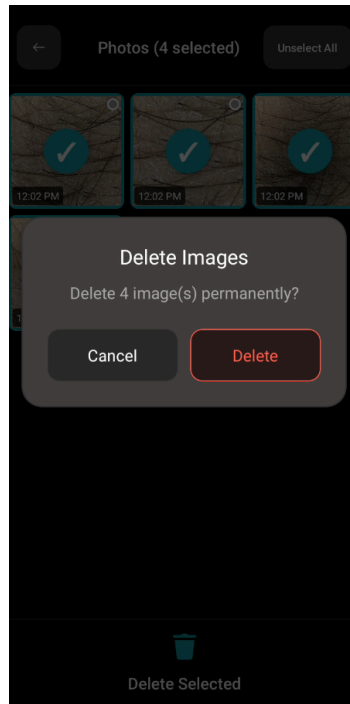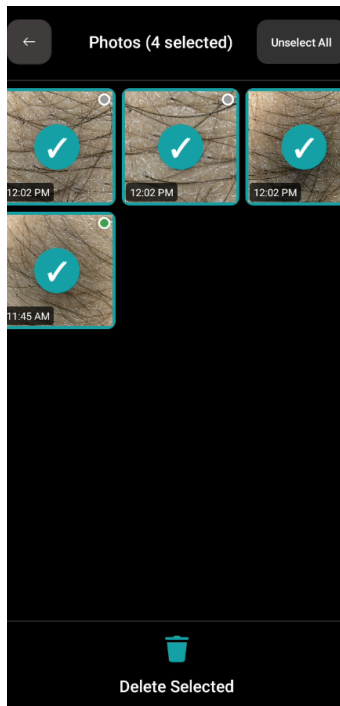


**Image Viewing:**

- View images on device screen.
- Access stored images via patient/case manager.
- Pinch to zoom on image for detailed examination.
- Zoom and review captured data.




#### 5.2.4. User settings:

- The user can set up Wi-Fi and date/time.
- The user can also log out via the logout button.

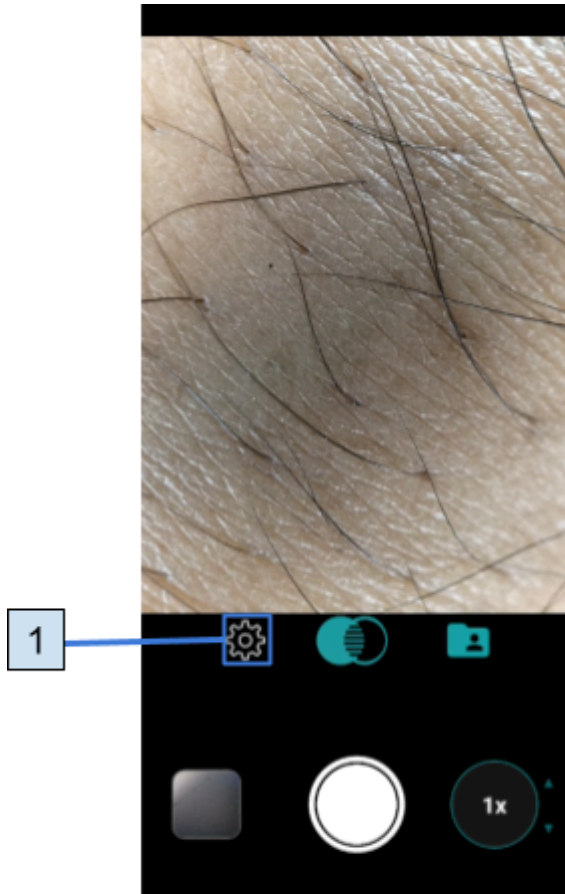

Camera screen

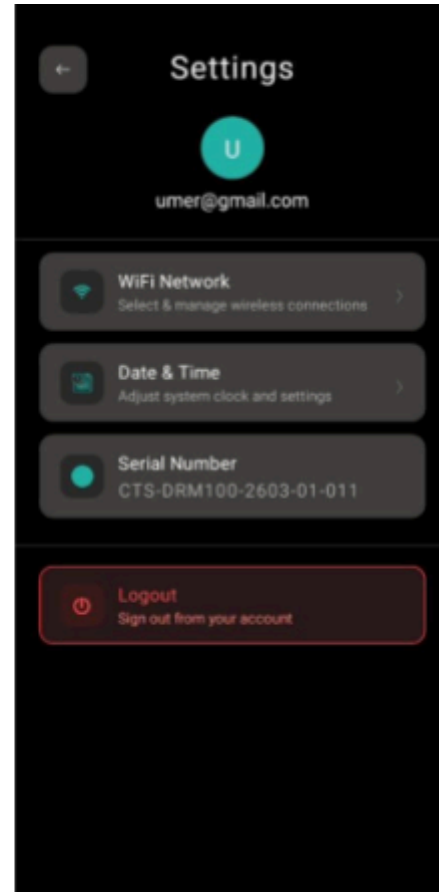


Settings

#### 1. User Settings


### 5.3. Wi-Fi Connectivity:

Automatically connects to previously saved networks.

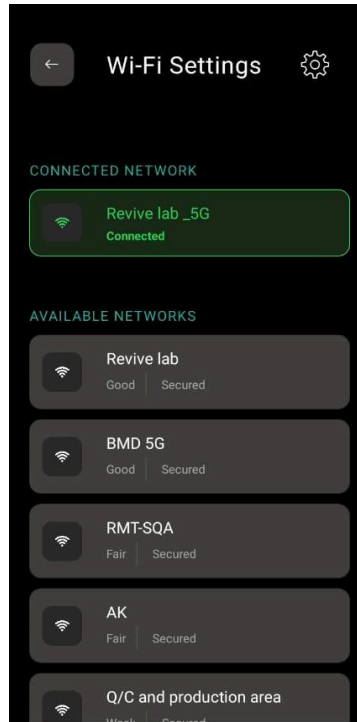


## 5.4. System Calibration (Date & Time):

- Accurate clinical documentation relies on precise chronometric data.
- The System Menu provides access to clock calibration.
- The Cutiscope Dermatoscope maintains internal chronometric consistency across power Cycles.

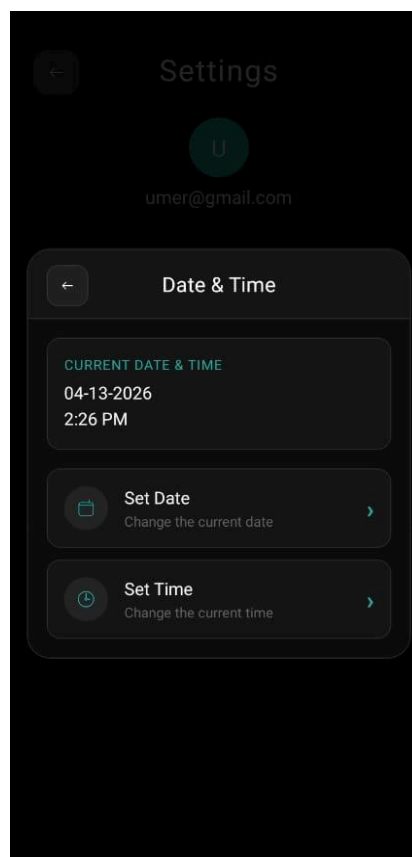

## 5.5. Power ON/OFF:


- When the device is Powered ON, upon pressing the power button, the power menu appears so the user can:
  - Power OFF
  - Restart
  - Lock/Standby
- If the device is already Powered Off then long press the power button to start the device.

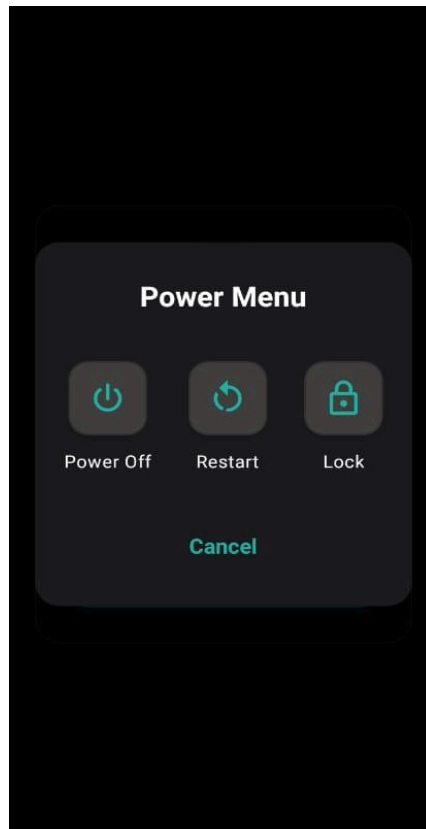



## 5.6. Sleep Mode:

The standby mode for the device is 2 minutes for no-touch input.

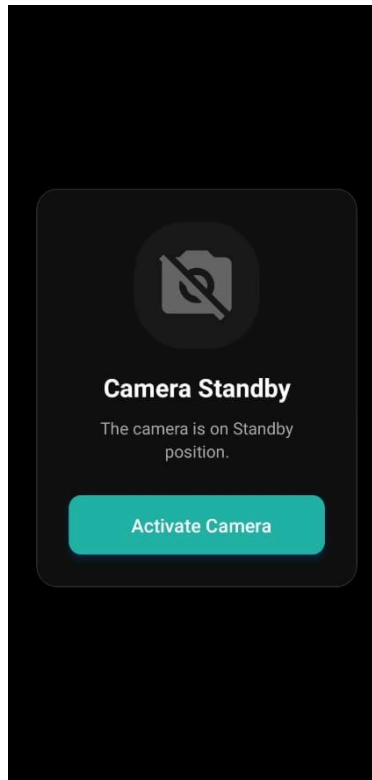

## 5.7. Battery Charging:


- Connect the USB cable to the device
- Charge via adapter
- The indicator shows the charging status.
- To ensure device longevity and data integrity, illumination is automatically switched off if the internal power level reaches a critical threshold (20%).

## 5.9 Web Dashboard:

The web dashboard is a secure, browser-based interface designed for clinicians to manage patient data and medical imagery. It connects to cloud storage and provides a central hub for viewing patients records, photos, and clinician activity.


- View the total number of patients, photos captured, and new patients added weekly.
- Access organized folders named by patient ID or name, and easily find patients using search by name or ID, with details such as creation date and photo count displayed.
- See login history and user actions for security and auditing.
- Each clinician uses an email and password to access or create their account.

### Operating Instructions:

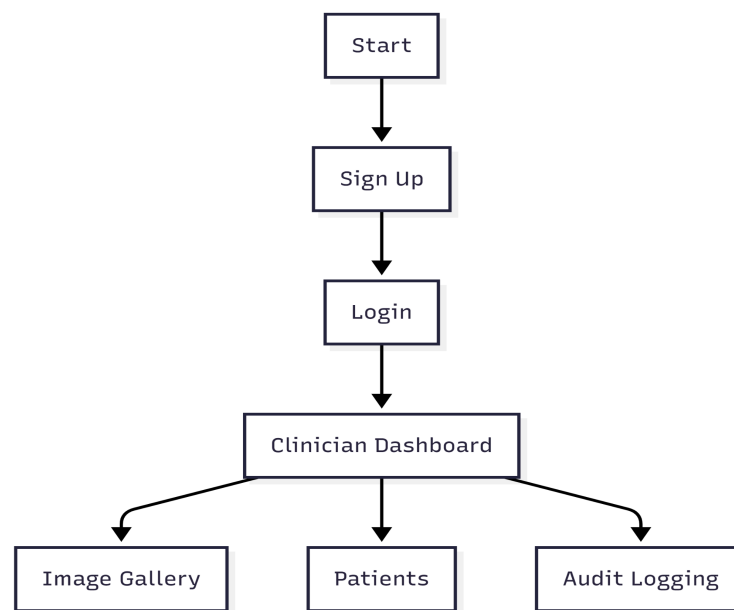


#### 5.9.1. Sign Up /Log In:

- Visit the site: 35.154.32.201:4040/login
- The user on first use will sign up and then using the sign-up credentials he/she can login whenever they want to.

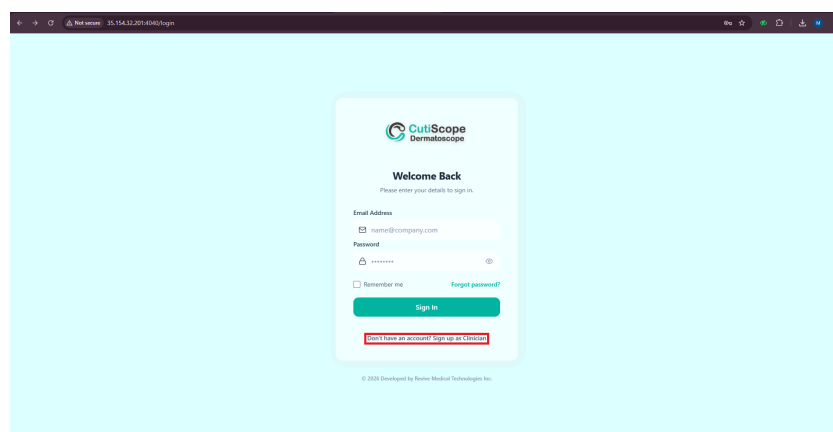


## 5.9.2. Dashboard:

After sign-up/login, the dashboard appears:

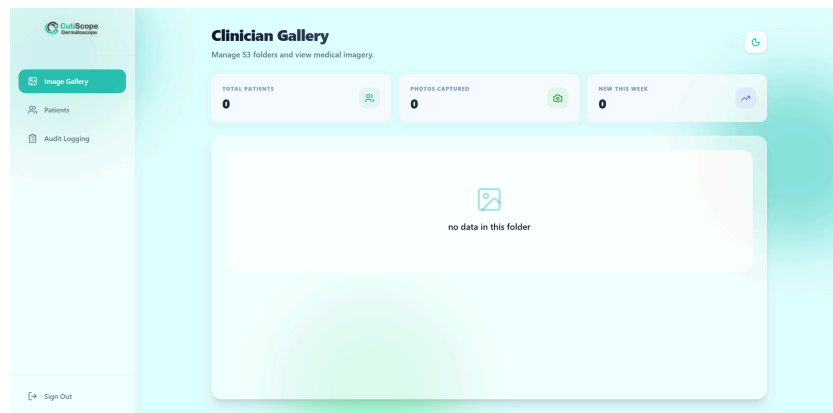


### 5.9.2.1. Image Gallery:

- In the image gallery, the Clinician can see all the images that he/she captured with the device.

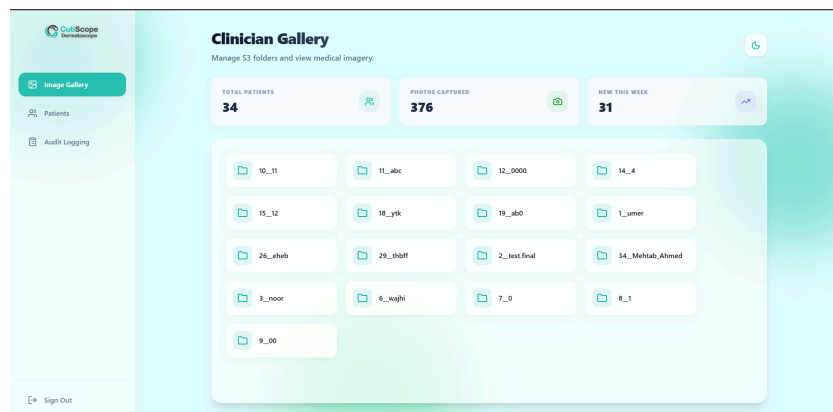


### 5.9.2.2. Patients:

1. In the patients tab, the clinician can add the new patient.
2. The list of patients that are added can also be seen.

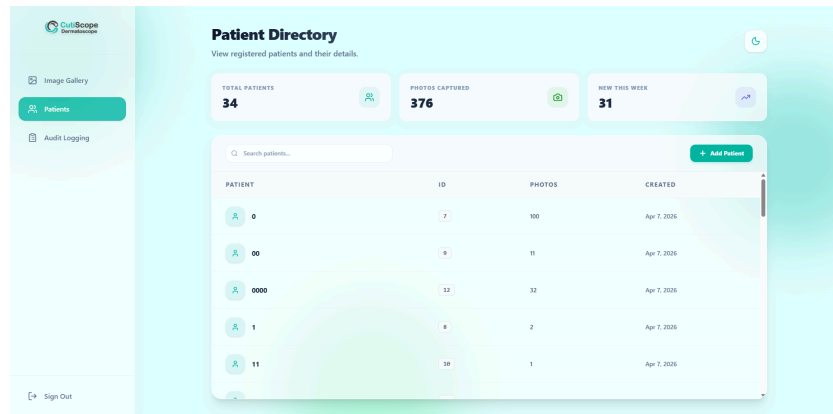

### 5.9.2.3. Audit Logging:

- Every time you log in, whether from the web dashboard or a device, the date and time are recorded.
- Each time you add a new patient, the system saves the patient's name either from device or web dashboard along with the exact date and time of the entry.

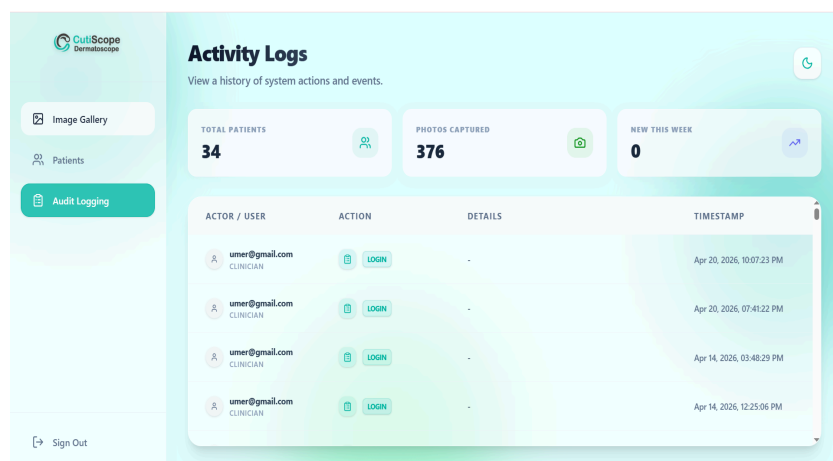

## 6. Electrical Safety:

- Do not attempt to open, modify, or repair the device, as this may expose internal electrical components.
- Use only the specified power adapter and the charging cable provided or recommended by the manufacturer.
- If the device shows signs of battery swelling, leakage, or overheating, discontinue use immediately.
- Keep the device away from strong electromagnetic sources such as high-frequency equipment.

## 7. Cleaning Protocol (Routine - Before/After Use)

| Area                               | Frequency          | Procedure                                                                                                                                                                                                                    |
|------------------------------------|--------------------|------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| <b>Lens      Contact<br/>Plate</b> | Before & After Use | <ul style="list-style-type: none"> <li>• Wipe with a lint-free cloth using non-alcohol wipes on external surfaces; allow to air dry.</li> <li>• Gently clean with a soft optical cloth; avoid excessive pressure.</li> </ul> |
| <b>Device      Outer<br/>Body</b>  | Daily / As Needed  | <ul style="list-style-type: none"> <li>• Clean with a slightly damp cloth and disinfect with Non Alcohol Disinfectant.</li> </ul>                                                                                            |
| <b>Buttons</b>                     | Daily              | <ul style="list-style-type: none"> <li>• Wipe gently to remove dust and contaminants.</li> </ul>                                                                                                                             |
| <b>Charging Port</b>               | Weekly / As Needed | <ul style="list-style-type: none"> <li>• Inspect and clean using dry air or a soft brush (no liquid).</li> </ul>                                                                                                             |
| <b>Full      Device<br/>Check</b>  | Daily              | <ul style="list-style-type: none"> <li>• Ensure no residue, moisture, or visible contamination remains.</li> </ul>                                                                                                           |

## 8. Troubleshooting Guide

| Problem                        | Possible Cause           | Check / Corrective Action                                                                                                       |
|--------------------------------|--------------------------|---------------------------------------------------------------------------------------------------------------------------------|
| <b>Device not turning ON</b>   | Battery drained          | <ul style="list-style-type: none"> <li>Recharge device using a USB-C charger</li> </ul>                                         |
|                                | Internal fault           | <ul style="list-style-type: none"> <li>Contact manufacture (refer to <a href="#">Section 9 Manufacturer Details</a>)</li> </ul> |
| <b>Poor image quality</b>      | Dirty contact plate      | <ul style="list-style-type: none"> <li>Clean using proper cleaning procedure</li> </ul>                                         |
|                                | Improper positioning     | <ul style="list-style-type: none"> <li>Adjust focus and placement</li> </ul>                                                    |
|                                | Wrong mode selection     | <ul style="list-style-type: none"> <li>Switch illumination mode</li> </ul>                                                      |
| <b>Device not charging</b>     | Loose connection         | <ul style="list-style-type: none"> <li>Check adapter, cable and port</li> </ul>                                                 |
|                                | Damaged adapter/cable    | <ul style="list-style-type: none"> <li>Replace or service</li> </ul>                                                            |
| <b>Overheating</b>             | Prolonged use            | <ul style="list-style-type: none"> <li>Turn OFF and allow cooling</li> </ul>                                                    |
|                                | High ambient temperature | <ul style="list-style-type: none"> <li>Use in recommended environment</li> </ul>                                                |
| <b>Software not responding</b> | Software glitch          | <ul style="list-style-type: none"> <li>Restart application</li> </ul>                                                           |
| <b>Image not saving</b>        | Storage full             | <ul style="list-style-type: none"> <li>Free up storage space</li> </ul>                                                         |

|                               |                                    |                                                                                                                                                       |
|-------------------------------|------------------------------------|-------------------------------------------------------------------------------------------------------------------------------------------------------|
| <b>WIFI Not Connecting</b>    | Weak signal                        | <ul style="list-style-type: none"> <li>• Restart the device</li> <li>• Move closer to the router</li> <li>• Check if networks are detected</li> </ul> |
| <b>Lens LED Not Working</b>   | Mode switching issue               | <ul style="list-style-type: none"> <li>• Switch polarization modes</li> </ul>                                                                         |
|                               | Low battery protection (below 20%) | <ul style="list-style-type: none"> <li>• Charge the device fully</li> </ul>                                                                           |
| <b>System Standby</b>         | Device is in power save mode       | <ul style="list-style-type: none"> <li>• Press the power button to turn ON the device.</li> </ul>                                                     |
| <b>Initialization Refresh</b> | Camera requires calibration        | <ul style="list-style-type: none"> <li>• Tap the “Wake Up” prompt.</li> </ul>                                                                         |



## 9. Manufacturer Details

**Manufacturer:** Revive Medical Technologies Inc.

**Address:** Building 2A, W1 Street, Rawat Industrial Estate, Islamabad, 46220

**Phone (Pakistan):** +92 3478807708

**Phone (USA):** +1 (707) 5618 771

**Email:** [info@rmt-usa.com](mailto:info@rmt-usa.com)

**Website:** <https://rmt-usa.com/>

**Country of Origin:** Pakistan

For technical support or service inquiries, please contact the manufacturer using the details provided above.